
Atlanta

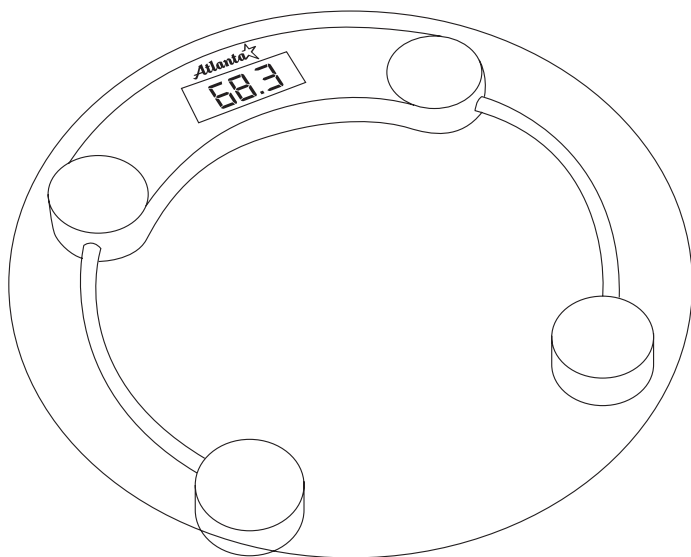ATH-6142

Personal electronic scales

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
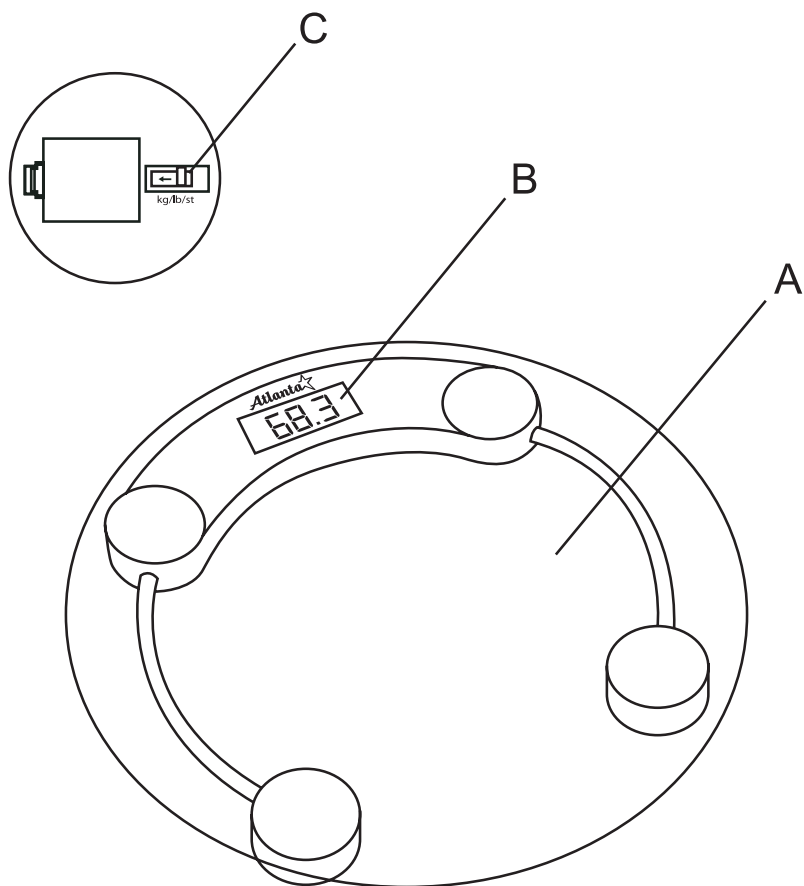
INSTRUCTION MANUAL

DEAR CUSTOMER!


CONGRATULATIONS ON YOUR GOOD CHOICE! NOW YOU ARE THE OWNER OF ATLANTA™ PRODUCT NOTABLE FOR ITS PROGRESSIVE DESIGN AND GOOD QUALITY PERFORMANCE. WE HOPE OUR PRODUCTS WILL STAY WITH YOU FOR YEARS ON END.

BEFORE USING

Please, attentively read this instruction before the beginning of operation with scales. Insert two AAA batteries as shown in the correct orientation when installed (batteries not included).

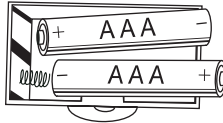If scales are delivered with the established CR2032 battery, it is enough to remove blocking film from the battery, as is shown in drawing to work with scales.

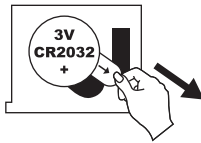For battery replacement, open a cover and establish a new battery observing thus polarity.

SAFETY MEASURES

- Electronic scales are designed for operation from the power supply with characteristics are specified on its marking.
- THIS DEVICE MAY NOT BE USED BY PERSONS (INCLUDING CHILDREN) WITH REDUCED PHYSICAL, SENSORY OR MENTAL CAPABILITIES OR LACK OF EXPERIENCE AND KNOWLEDGE IF THEY ARE NOT SUPERVISED OR HAVE NOT BEEN INSTRUCTED CONCERNING THE USE OF THE DEVICE BY A PERSON RESPONSIBLE FOR THEIR SAFETY.
- DON'T LEAVE THE DEVICE UNATTENDED IF IT IS SWITCHED ON OR WHEN CHILDREN ARE PRESENT NEAR IT. DON'T LET THE CHILDREN PLAY WITH OR USE THE DEVICE. CHILDREN AT THE AGE OF 8 AND OVER ARE ALLOWED TO USE THE DEVICE IF THEY ARE CONSTANTLY SUPERVISED BY ADULTS.
- In order to reduce the risk of damage of the device:
 - do not open up the housing of the product;
 - do not expose the electrical part of the device to rain or moisture;
 - do not immerse the device in water or other liquid;

EN
3-7RUS
8-12KZ
13-17

- The glass surface of scales can be very slippery if your feet damp or you use socks from a slippery material or footwear with a slippery sole. In order to avoid accidents, do not use scales in a bathroom and become on scales carefully.

Never use the device on a rough or slippery surface.

- The device is designed for operating indoors, for domestic use only. Commercial use of the device is prohibited.

If in an operating time there were any malfunctions, do not use scales. Address in the nearest service centre with the warranty coupon.

- Do not place the device close to any sources of heat, such as: heating radiators, gas stoves; or in areas subject to direct sunlight, excessive dust, mechanical vibration or impact loads.

- After transportation under low temperature conditions, the device can be connected to the mains only after two hours in a warm room.

- Do not use any chemically active liquids (petrol, kerosene, alcohol, various solvents) for cleaning the housing of the device.

- Any faults in storage, transportation or operation of the device are not covered by the warranty.

BASIC ELEMENTS OF SCALES

- A - case
- B - liquid crystal display (LCD)
- C - "Unit" switch (option)

PERFORMANCE SPECIFICATIONS

Max. capacity:	180 kg
Graduation:	0.1 kg
Power:	1x3 V lithium battery, type CR2032 or two AAA batteries
Service life:	5 years

With due care in the course of operation the mechanical life of the device may considerably exceed the mechanical life specified by the manufacturer.

Date of production is represented on the label of the product in the format MM.YYYY

SCALES OPERATION

1. Establish scales on a firm, steady, horizontal surface. If the scales are equipped with the UNIT switch you can set form weighing according to convenient for you.

2. Will rise on scales, the display switch on and scales start to measuring. Try not to move, your weight will be reflected in the display in kg. For example, 68 kg of 300 gr.
3. Scales are intended for weighing of weight to 180 kg. In case of overweight the indicator will display an error.

4. After the termination of weighing scales will be disconnected automatically in 10 seconds.

EN
3-7RUS
8-12KZ
13-17

NOTES

- Scales are supplied by the indicator of low batteries. If the inscription is reflected in the display «LO», you should replace a battery in scales.

- That weighing was more exact, it should be produced in the same place of a room. Try to stand thus directly.
- Do not store scales in vertical position with battery as in such position scales can automatically will join and to plant the battery.
- Do not leave a battery in scales if you are not going to use scales long term.
- If as a result of double weighing scales have shown different results, your true weight is between two these values.
- The manufacturer reserves the right to itself to make to a electronic scales design.

CLEANING AND LEAVING

- Scales can be wiped slightly damp rag.
 - Do not use chemically active or abrasive washing-up liquids for cleaning of scales as it can damage their surface.
 - Do not clean scales by means of tools.
 - Never wet, and do not immerse an item in water or other liquid.
- OFTEN MEETING MALFUNCTIONS, THEIR REASONS AND ELIMINATION METHODS

INSTRUCTIONS FOR ENVIRONMENTAL PROTECTION

The old appliances are the assemblage of technical materials and therefore cannot be utilized with household rubbish! For this reason, we would kindly ask you to support us in resource saving and environmental protection and deliver this appliance into a utilization reception center (if there is one).

Troubleshooting	Possible cause	To solve the problem
Scales do not work	Electronic or static error	Remove a battery, wait some seconds, and then establish it into place. Make sure that the battery was not discharged differently replace the battery
On the display the inscription «LO» is displayed	The battery is low	Replace a battery
On the display the inscription «Err» is displayed	The peak load of scales, component of 180 kg is exceeded	On scales it is possible to measure only that weight, which does not exceed the established peak load
On the display the inscription «Err» is displayed	Scales were not ready to measurements (on the display "0.0" have not been established)	To repeat measurement, having waited zeroings of indications on the display
On the display the inscription «C» is displayed	Weight read with error	To repeat measurement, having waited zeroings of indications on the display

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