# SCARLETT®

# SC-BS33ED110

## INSTRUCTION MANUAL

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**GB DESCRIPTION**


1. Body
2. Display
3. Electrodes
4. 2 batteries 1.5V AAA ( included)
5. Batteries compartment

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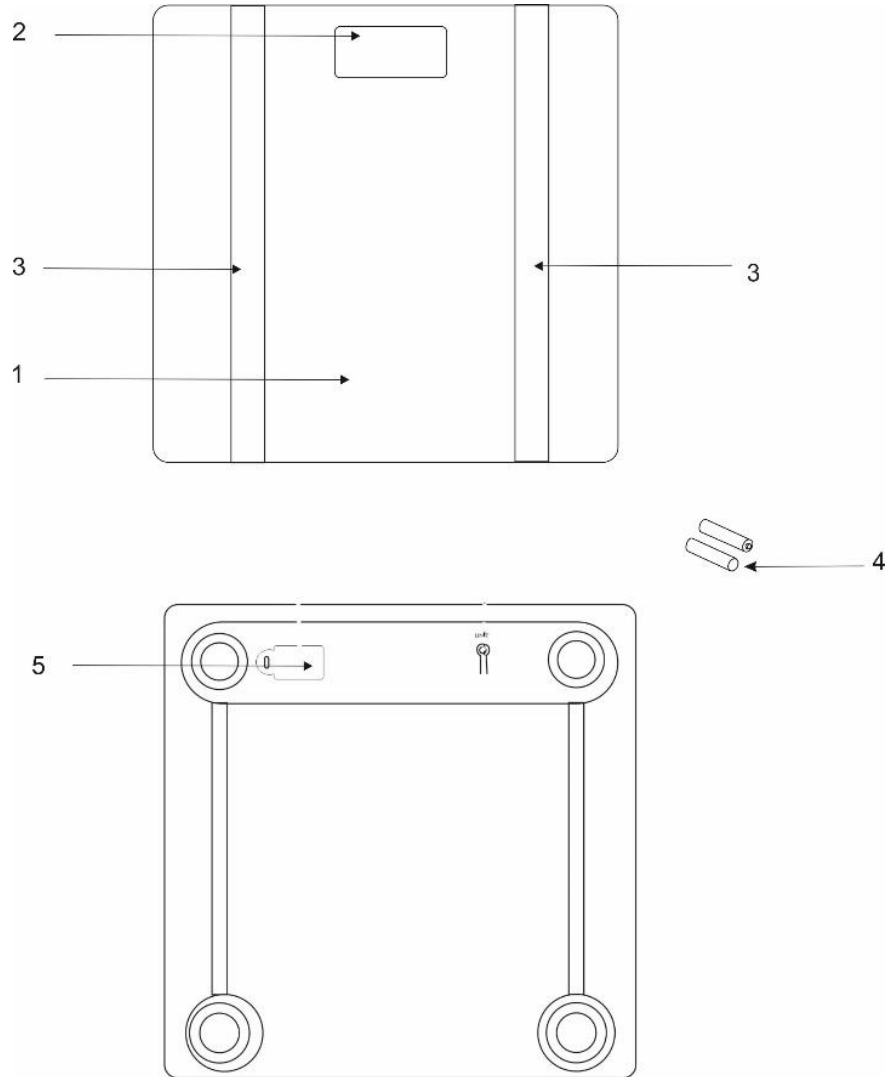## **GB INSTRUCTION MANUAL**

### **IMPORTANT SAFEGUARDS**

- To avoid breakdown during use read this guide carefully before operating the appliance. Mishandling can cause breakdown of the appliance.
- Only for domestic use. The appliance is not designed for industrial or commercial use.
- Do not expose the scale to shocking.
- Do not jump on the scale.
- Do not overload the scale.
- It is forbidden for use by persons with medical implants (for example, with pacemakers). Otherwise, faulty operation of implants might occur.
- Do not use during pregnancy.
- Attention! Do not step on the scale if your feet or the scale surface are wet – you might slip.
- Batteries may be life-threatening if swallowed.
- Store the batteries and the scale out of the reach of children. The person who swallowed a battery, should immediately consult a doctor.
- Do not allow wrapping material getting into the hands of children (choking hazard).
- It is forbidden to recharge or restore batteries in other manner, disassemble them, or throw into fire.
- Measuring body fat percentage and other parameters in children below 10 years of age may yield invalid results.
- Diagnostic scale may be used for simple weighing without determining water, body fat or muscle mass parameters in children of any age. Minimum weighing limit is 2 kg.
- Unreliable results of measuring body water, body fat or muscle mass might occur in professional athletes, persons prone to oedema, persons with irregular leg-to-body ratio.
- Never try to repair this appliance by yourself. In case of failure, contact the nearest seller.
- If the product has been exposed to temperatures below 0°C for some time it should be kept at room temperature for at least 2 hours before turning it on.
- The manufacturer shall reserve the right to introduce minor changes into the design of the scale, which do not affect its safety, performance capacity or functionality of the appliance.
- Manufacturing date is indicated on the unit and/or on packing as well as in accompanying documents in the following format XX.XXXX where first two figures XX is the month of production, and following four figures XXXX is the year of production.

### **BATTERIES USAGE RULES**

- Battery contains harmful substances for human body. Keep the wrapping materials out of reach of children (choking hazard). The person who swallowed a battery, should immediately consult a doctor.
- Change the batteries regularly. Change all batteries at once using the elements of the same type.
- It is forbidden to recharge or restore batteries in other manner, disassemble them, or throw into fire or hot-wire.
- Discharged battery can be a reason of scale breakdown. Remove the batteries compartment in case of long out-of-action time.
- If one of the batteries is discharged, put protection gloves and clean the batteries compartment with dry cloth.
- Batteries may contain toxic substances harmful for health and environment. The batteries should be disposed in strict compliance with existing law. Do not dispose the used batteries to the waste bin.
- Do not throw the batteries into fire! Highly explosive!

### **PRINCIPLE OF OPERATION**

- The principle behind this scale is bioelectrical impedance analysis (B.I.A.). The measurement of the percentage of different tissues occurs within seconds with electric current, which passes through high-effective sensor on scales platform and which is not perceptible for a human, safe and harmless.
- It is important that during measurement you take into account the following: measurement of the body fat percentage must be performed only barefoot; it is advisable that you slightly moisturize the soles of your feet. Absolutely dry soles might be the cause of an unsatisfactory result, as they possess extremely low electric conductivity.
- If, when measuring the percentage of fat, muscle tissue, bone and water content, the obtained sum of all components has become greater than 100%, this is not a sign of scale defect. It can be explained from the physiological point of view. Operation of the diagnostic scale is based on the analysis of bioimpedance, which determines the amount of fat and water in muscle and bone tissues themselves. Thus the variation of different parameters should be taken into account separately from each other, without summarizing the obtained results.
- This model is working with Bluetooth, compatible with smartphones on IOS and Android platforms. Use "HOME INTELLECT" program for convenient checking and fixation of measurements.

### **HELPFUL HINTS**

- Make sure the scale places on a flat, even and firm level surface. Do not put the scale onto a carpet or soft surfaces.
- To ensure correct weight measurement it is advisable that the scale stand in the same place. By moving the scale you risk putting it on an uneven surface, which might affect the correctness of your weight reading.
- To ensure correct weight measurement you need to weigh at the same time, without clothes or footwear, and before eating.
- To ensure precise weight measurement it is not recommended to weigh earlier than 2 hours after waking.
- Step on the scale. Do not lean against anything, stand straight.

#### LIMITATIONS

- The measuring of fat percentage and other measurements may be unreliable and inexact in the following cases:
  - For children younger than 10 years old
  - for professional sportsmen and bodybuilders;
  - for pregnant women;
  - for people with hyperthermia, at dialysis treatment, edema, osteoporosis;
  - for people taking cardiovascular pills (with cardio and vascular abnormalities)
  - for people taking vasohypotonic and vasohypertonic pills

#### GETTING STARTED

- Insert type AAA batteries observing polarity into a special compartment at the bottom of the scale. Close the compartment. The scale is ready to operate.

#### WEIGHING

Note! When the scales are turned on, the ambient temperature is displayed in the lower right corner of the display

- Input your personal data to get a precise measuring body fat percentage and other body measurements
- The scale has a memory for 10 persons data which can be used by your family.
- Connect the Internet using your smartphone(WI-FI, 3G or 4G are recommended)
- Switch on Bluetooth function.

#### CONNECTING THE DEVICE TO THE HOME INTELLECT APP


<https://home-intellect.ru/appdownload/>

#### 1. Install the HOME INTELLECT mobile application

- 1.1. Go to the Home Intellect website using the links above.
- 1.2. Install the Home Intellect app on your mobile device or tablet.
- 1.3. Make sure your smartphone or tablet is connected to 4G internet or Wi-Fi.
- 1.4. Make sure your smartphone or tablet is running iOS 8.0 (or higher) or Android 4.4 (or higher).


#### 2. Create an account in the mobile application:

- 2.1. Open the application and create a new account in it:  
Enter your email address

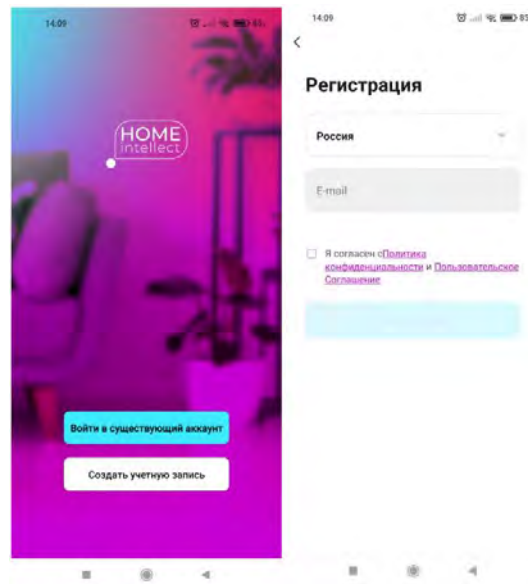
An automatic email with a confirmation code will be sent to the entered email address, enter this code in the field that appears

- 2.2. Set the desired password for your account
- 2.3. Your new account will consist of an email and a password you created
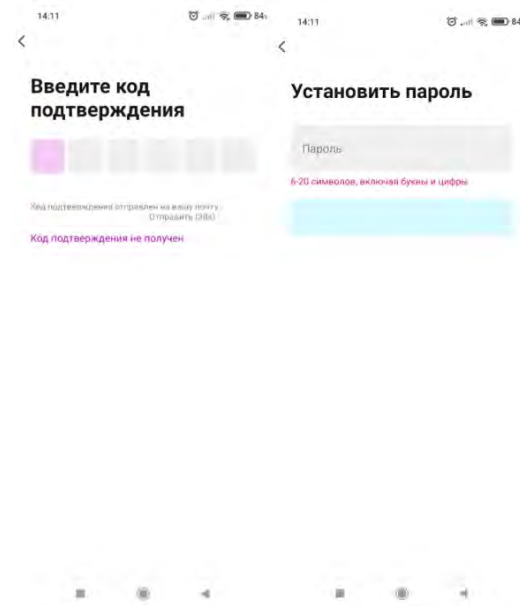
### 3. Connect the device to the mobile application


- 3.1. Make sure geolocation is enabled on your device.
- 3.2. Click on the plus icon in the top corner of the screen.
- 3.3. Select "Add Device"
- 3.4. Next, you need to turn on the Bluetooth function on your smartphone - the application will automatically detect the scales (we recommend standing on the scales to turn them on)
- 3.5. Select «Добавить»

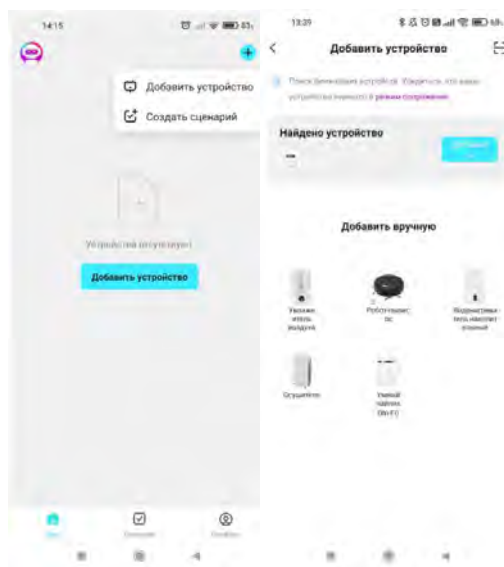- 3.6. When you go to the next screen, click on the plus icon - the application will take no more than 2 minutes to add the device.

## Main features of the application:

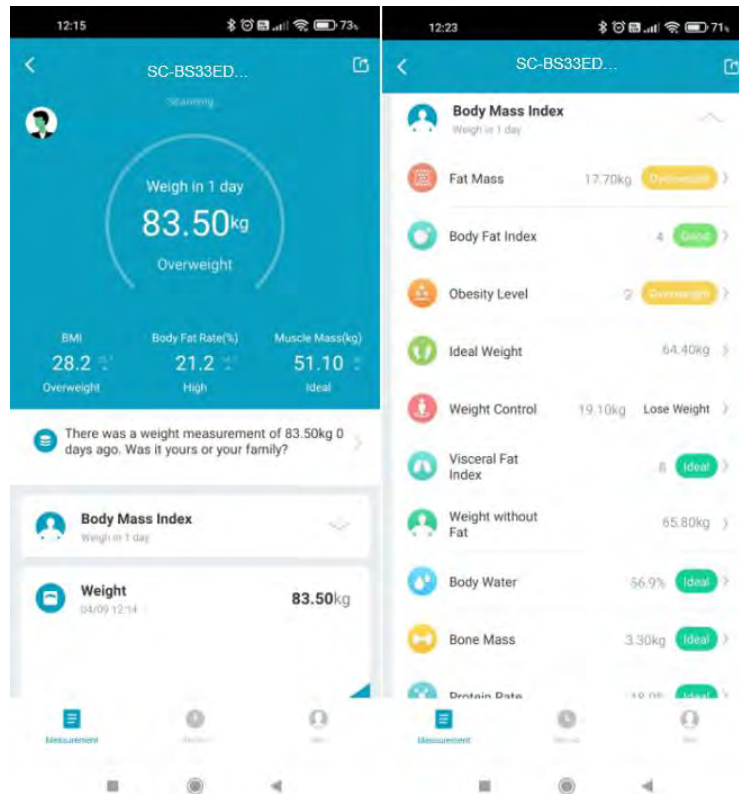- Measurement of total body weight
- Body mass index measurement
- Measuring the level of fat mass in kilograms
- Body fat mass index measurement
- Assessment of overweight
- Receive ideal weight recommendations
- Receive weight management advice
- Measurement of visceral fat mass index
- Calculation of lean mass
- Indication of the level of water in the body
- Measurement of bone mass in kilograms
- Protein measurement
- Measurement of basal metabolic rate in kilocalories
- Estimating your body age
- Description of your body type
- Calculation of the overall score of physical indicators
- Maintaining the history of weighing

### WHAT TO DO IF THE APP DOES NOT DETECT THE DEVICE:

- 1) Make sure your mobile device has Bluetooth enabled
- 2) Press and hold the change unit button on the bottom of the device until DEL appears on the screen - after that the device will reboot to reconnect to the mobile application
- 3) Make sure this unit is not paired with another mobile device. In this case, remove it from the application on another mobile device.

*The appearance of the mobile application, its interface, icons and the location of functions may change without prior notice to improve its performance.*

### AUTOMATIC RESET FUNCTION

- Your scale has an automatic reset function. After the data is displayed on the screen, the data will be reset in a few seconds and the screen will show 0.0 again. You may conduct the weighing again or start the weighing by setting a new user.

### CLEANING AND MAINTENANCE

- Wipe the scale with a soft cloth using a detergent and wipe dry.
- Do not use organic solvents, aggressive chemicals or abrasives.

## STORAGE

- Observe all requirements from the CLEANING AND MAINTENANCE Section.
- Make sure there are no objects on the scale during storage.
- Store the scale in a cool dry place.

The symbol on the unit, packing materials and/or documentations means used electrical and electronic units should not be tossed in the garbage with ordinary household garbage. These units should be passed to a special receiving point.

The symbol on the unit, packing materials and/or documentations means the battery used in this product should not be tossed in the garbage with ordinary household garbage.

For additional information about the actual system of the garbage collection, address to the local authority. Valid utilization will help to save valuable resources and avoid negative work on the public health and environment which happens with incorrect using garbage.

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